



slee solutions
for
pathology


MSM

CAROUSEL SLIDE STAINER

DESIGN &
MANUFACTURING
MADE IN GERMANY

MSM Carousel slide stainer

With its simple user interface and sophisticated technology, the **MSM** carousel stainer is designed to meet the requirements for routine staining in histology and cytology laboratories. The device is equipped with 24 incubation stations containing up to 400 ml. The flex-programming option allows different staining protocols while using the same reagents. An integrated UPS (uninterruptible power supply, with built-in batteries) ensures highest process reliability even in environments with frequent power failures. A reliable extraction system with carbon filter reduces contact with reagent vapors for increased workplace safety.

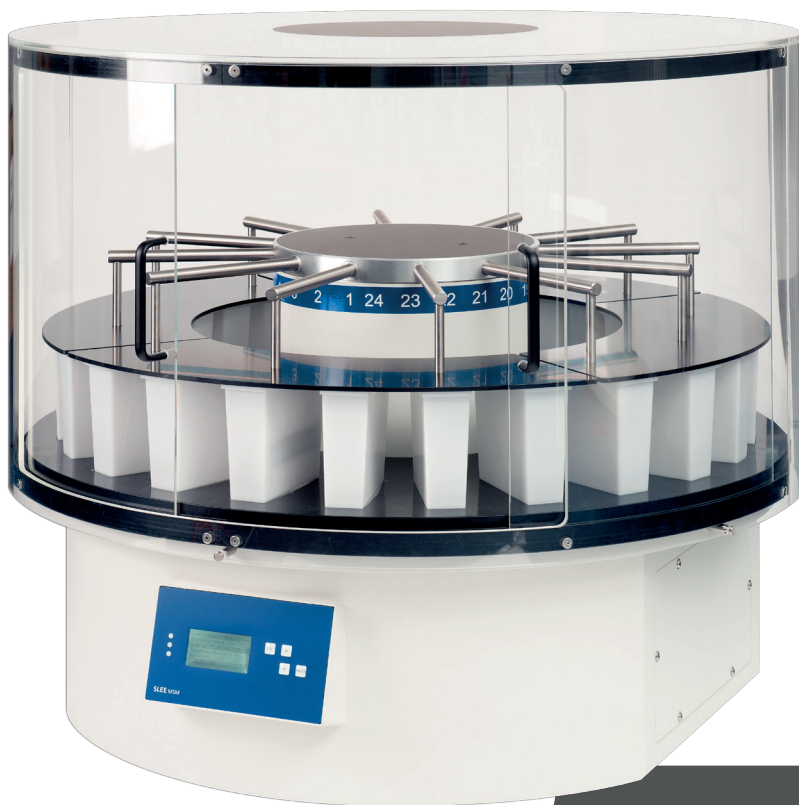

Highlights

- Easy loading and full process control due to rotating table principle
- Integrated UPS system protects your samples
- Flex-programming
- System status monitor
- Password protection
- Intuitive control panel

Item No.**MSM 100 - 240 V, 50/60 Hz**

12000205

Standard components

Staining cuvettes (plastic, 24 pcs.)**Transport basket** (plastic, 1 pcs.)
capacity 30 standard size micro slides**Mains cable****Manual**

Specifications

Staining cuvettesmaterial plastic
dimensions (W x D x H) 130 x 40 x 135 mm
volume 400 ml**Transport basket**material plastic
capacity 30 standard size micro slides**Programming** 20 staining programs,
freely programmable**Electrical parameters** 100 - 230 V, 50/60 Hz**Dimensions** (W x D x H) 870 x 870 x 750 mm**Weight**

w/o accessories 80 kg

Options

Item No.**Hood** for MTP / MSM (Perspex) 11000220**Staining cuvettes** 12000260**Slide basket** (30 slides) 12000261**Drying station** 12000256**Exhaust fan** 11000210**Exhaust fan**
incl. carbon filter 11000211**Carbon filter** 11000254



SLEE medical GmbH · Am Neuberg 14
55268 Nieder-Olm · Germany
www.slee.de

T: +49 (0) 6136 76997-0
E: mail@slee.de



www.slee.de



YouTube



LinkedIn



Instagram

